**ACTIVE
LUXE**

H · A · I · R

LED
Light Therapy
Cap

User Manual
Uporabniški priročnik
Korisnički priručnik


1

2

Working Principle

Before getting started, please take a look at this user manual to get tips for usage and maintenance. This product delivers medical-grade irradiance via red light and infrared light. The 660 nm and 850 nm (invisible) wavelengths have been proven to be effective in pain relief, wound healing, inflammation reduction, etc. Consistency is the key. For the best results, we recommend using 20 minutes per time, 3-4 times per week. If you're new to RLT, please start with a 2-minute treatment, and slowly work up to a 10-minute session for 2 weeks as your body adapts to the therapy. Results can be seen after 8-12 weeks of consistent use. 660 nm Red Light could speed up the metabolic reaction of cells, promote skin health, and speed up healing. Help the hair regrowth.

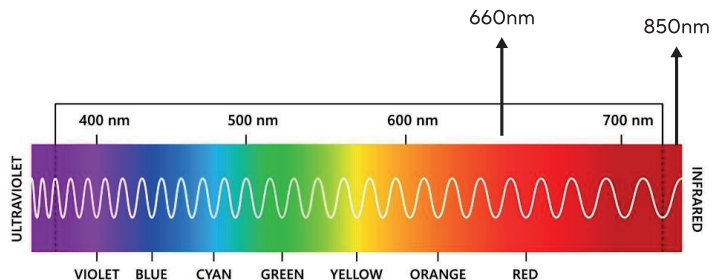850 nm Near Infrared light could penetrate deeper into the tissue, organs, and joints, making it ideal for enhancing muscle recovery and reducing joint pain. Helps the hair regrowth. Note: 850 nm is invisible to the human naked eye.

Technical Specifications

Model ID	T100Pro
Working Voltage	DC 5 V
Power	10 W
LED Source	5050SMDLED
LED quantity	150 PCS
Wavelength	660:850=1:2
Dimension	30 × 19 × 16 cm
Timer	30 minutes AutoTurn OFF
Diming	Level-1 to Level-3
Battery Capacity	2000 mAh

Packaging Contents

1. Cap packing box
2. Red Light Therapy cap
3. USB charging wire
4. User Manual

Using Steps

1. Use the charger to fully charge the cap.

4. Adjust the cap size by adjusting the cap buckle.

2. Put the inner layer in the cap.

5. Wear the cap on your head.


3. Fasten the Velcro along the edge of the cap.

6. Fasten the Velcro along the edge of the cap.

Operating Instruction

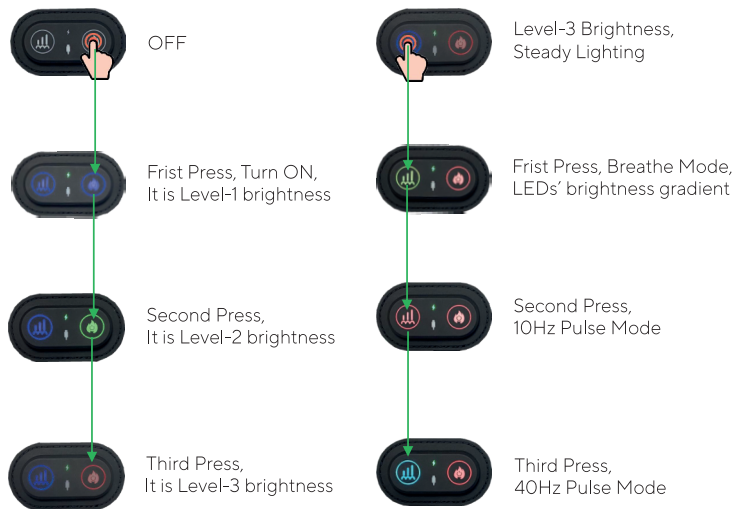Default Timer:
30 Minutes

Long press on this key to turn off the Wrap.

Warning

1. The product should be used while sitting, do not walk while using the product.
2. Do not use this device on infants.
3. The device is not waterproof, please keep it away from liquid or essential oils.
4. Only use the provided power cable.
5. Consult your doctor prior to use if you are: pregnant; photosensitive; recently injected steroids or use topical steroids.
6. Do not use if you are allergic to neoprene.
7. Do not use liniments, salves or ointments that contain heat producing ingredients as skin burns or irritation may occur.
8. Do not use this device if you are sensitive to light, or you are currently taking medications that increase light sensitivity. Doing so may result in a photosensitive reaction. Discontinue use immediately if you experience any discomfort or sensitivity.
9. Do not use on open wounds, heat may increase circulation and cause bleeding, and warmth may increase infection in some climates.
10. Contact your doctor before purchasing or using this device if you have any concerns about using this device in general or while taking any current medications.

Maintenance

Make sure to turn off the device and unplug the power adapter. Dot fold the wrap, keep it flat when you not using it.

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